

STOMICS FILE FORMAT MANUAL

Revision History

Manual Version: A0
Software Version: V4.1.0
Date: Apr. 2022
Description: Initial release

Note: Please download the latest version of the manual and use it with the software specific for this manual.

©2022 Beijing Genomics Institute at Shenzhen (BGI-Research). All rights reserved.

1. The products shall be for research use only, not for use in diagnostic procedures.
2. The Content on this manual may be protected in whole or in part by applicable intellectual property laws. BGI-Research and / or corresponding right subjects own their intellectual property rights according to law, including but not limited to trademark rights, copyrights, etc.
3. BGI-Research do not grant or imply the right or license to use any copyrighted content or trademark (registered or unregistered) of us or any third party. Without our written consent, no one shall use, modify, copy, publicly disseminate, change, distribute, or publish the program or Content of this manual without authorization, and shall not use the design or use the design skills to use or take possession of the trademarks, the logo or other proprietary information (including images, text, web design or form) of us or our affiliates.
4. Nothing contained herein is intended to or shall be construed as any warranty, expression or implication of the performance of any products listed or described herein. Any and all warranties applicable to any products listed herein are set forth in the applicable terms and conditions of sale accompanying the purchase of such product. BGI-Research, Shenzhen makes no warranty and hereby disclaims any and all warranties as to the use of any third-party products or protocols described herein.

TABLE OF CONTENTS

CHAPTER 1: OVERVIEW

1.1. About Software	1
1.2. About Manual	1
1.3. Terminologies and Concepts	1

CHAPTER 2: FILE FORMAT

2.1. BAM	3
2.2. Mapped CID List with Reads Count File	3
2.3. Gene Expression File	4
2.4. Gene Expression Matrix	7
2.5. Image Pyramid	8

REFERENCES	9
------------	---

CONTACT US	10
------------	----

CHAPTER 1

OVERVIEW

1.1. About Software

STOmics Analysis Workflow¹ (SAW) software suite is a set of pipelines that are bundled to position sequenced reads to their spatial location on the tissue section, and quantify spatial gene expression.

SAW download (Docker Hub): <https://hub.docker.com/r/stomics/saw>

SAW Github: <https://github.com/BGIResearch/SAW>

1.2. About Manual

This manual includes descriptions of key files format generated from SAW, which help users better understand and make use of information from analysis results.

1.3. Terminologies and Concepts

Table 1-1 Terminologies and Concepts

Abbreviation	Full Name	Description
SN	Serial Number	Unique ID for STOmics Gene Expression Chip.
RIN	RNA Integrity Value	RNA integrity value measures the RNA degradation degree to indicate the integrity of RNA and evaluate the quality of the RNA sample. RIN values range from 1 (totally degraded) to 10 (intact). In STOmics analysis, only tissue sample with a pre-measured RIN value greater than 7 should be used for further sequencing and bioinformatics analysis.
CID	Coordinate ID	Spatial position identifier, the artificially synthesized barcode sequence unique to each spot on the STOmics Gene Expression Chip.
MID	Molecular ID	Molecular identifier (same as UMI), the artificially synthesized sequence unique to each mRNA molecule captured from the sample which helps to differentiate the number of reads contributed by mRNA expression level due to amplification. Two copies of native transcripts from the same molecule captured on one DNB will result in two independent reads with the same CID but different MID. In contrast, two reads with identical CID and MID were originated from the same transcript but got amplified.
DNB	DNA Nanoball	DNA nanoball is the product of rolling-circle amplification (RCA) that is linearly amplified from the original circular single-stranded DNA template. DNB is the smallest capture unit on the STOmics Gene Expression Chip.
Bin	Bin	Bin (or Square Bin) is the analysis unit on the gene expression heat map. A bin is a fixed-sized square in which the expression value in this square is accumulated. Bins are not overlapped. The value followed by “Bin” represents the side length of the square. For bin 1, each DNB on the STOmics Gene Expression Chip is shown as a spot, which means one spot only contains the data from one DNB. Bin N means one spot on the heat map is an aggregation of data from N×N neighbor DNBs. For example, a spot of bin 100 covers data from 10,000 DNBs.

CHAPTER 2

FILE FORMAT

2.1 BAM

The BAM² file format is a binary format for saving sequence alignment and gene annotation data. SAW **mapping** BAM adds custom tags in the BAM optional field to record reads coordinates, CID and MID information. **count** BAM adds annotation information in the tag field. Custom tags are described in Table 2-1.

Table 2-1 BAM custom tags

Tag	Description
Cx:i	x coordinate of CID.
Cy:i	y coordinate of CID.
UR:Z	The hexadecimal representation of uncorrected binary-encoded MID.
XF:Z	Mapping region on the reference genome. Valid value: 0=EXONIC, 1=INTRONIC, 2=INTERGENIC.
GE:Z	Annotated gene name.
GS:Z	'+' or '-', indicating forward/reverse strand respectively.
UB:Z	The hexadecimal representation of count corrected binary-encoded MID.

Example of mapping BAM:

```

E100026571L1C009R00301275185      16      1      3000095 255      26M121066N74M
*      0      0      GGCTTTTTTTTTTTTTTTTTTTTTTTTTTTTTCTAAATATTGGGTTTTATTAGC
ACCATGATAACTGTATATTAATTTGCACTGACTGTCATAACAAAATAC      G+:GFFGGGFFGFFGFFGFFGFFG
FFFFCFGFCFGGGFFGFFGFFGFFGFFGFFGFFGFFGFFGFFGFFGFFGFFGFFGFFGFFGFFGFFGFFGFFGFFG
NH:i:1  HI:i:1  AS:i:88 nM:i:0  Cx:i:4826      Cy:i:11598      UR:Z:6FA29

```

Example of count BAM:

```

E100026571L1C002R00703943265      1040      1      3082766 255      11M132671N89M
*      0      0      CTGCTGCAGCTTTTTTTCTTTGAGATTTATTTTATGCTATGTGTATGGGT
ATTTTGCCTGCATATATGTCTATGCACCATGTGTGTGAGTGCTTGAG      FFFFECCGFCGDCGDFGDFEE@EEG
IBFGGCGFFGACGFCGFFDGDGFFFFFEGCDFCGFFGG@FFF=EFFDGGGGGFDGFFFGGGGFFGGGFFGGGDFG
NH:i:1  HI:i:1  AS:i:88 nM:i:0  Cx:i:7767      Cy:i:18052      UR:Z:7AE49
XF:i:0  GE:Z:Xkr4      GS:Z:-  UB:Z:79E49

```

2.2 Mapped CID List with Reads Count File

mapping pipeline outputs mapped CID list file with reads counts for each CID. This file stores CID coordinates and reads count for each coordinate. The list does not have header. The three columns are x coordinate, y coordinate and MID count.

Example of mapped CID list with reads count file:

```

14195      16619      1
19945      14450      2
14548      9438      1

```

2.3 Gene Expression File

Gene expression file (GEF) is a type of gene expression distribution visualization file particularly designed for Stereo-seq. GEF is a hierarchically structured data model that stores multiple gene expression matrices in different bin sizes.

Each GEF container organizes a collection of spatial gene expression matrices. It includes two primary data objects, Group and Dataset. A Dataset is a multidimensional array of data elements. Group object is analogous to file system directory which organizes Datasets and other Groups in hierarchies.

The first level of GEF includes three Group object: “geneExp”(required), “wholeExp”(optional), and “stat”(optional). Group “geneExp” contains Groups of gene spatial expression data in one or multiple bin size. Group “wholeExp” contains Datasets that record expression level and gene type count of each coordinate in one or multiple bin sizes. Group “stat” saves gene names, total MID count and spatial pattern enrichment score of each gene. “Attributes” in each Dataset records the key metrics of that Dataset. Check <https://www.processon.com/view/link/610cc49c7d9c087bbd1ab7ab#map> to get the schematic diagram of GEF. The field names and field data types are described in Table 2-2. SAW outputs three GEF files in the whole process, please check the Table 2-3 to find the description.

Table 2-2 Gene Expression FileText Fields Description

Attributes			
File Attributes	Data Type	Example	Description
version	uint32	1	Gene expression file format version.
geftool_ver	uint32	0,6,2	geftool version. It can be used as an individual tool to manipulate GEF files.
geneExp/binN/expression: Dataset “expression” is a 1D array which stores coordinates and MID counts of each gene in the bin size of N, aggregated by gene name.			
Dataset Attributes	DataType	Example (bin1)	Description
minX	uint32	59820	Minimum x coordinate in bin N.
minY	uint32	102086	Minimum y coordinate in bin N.
maxX	uint32	73040	Maximum x coordinate in bin N.
maxY	uint32	120539	Maximum y coordinate in bin N.
maxExp	uint8/uint16/ uint32	28	Maximum MID count in a spot when the bin size is N. Data type for “maxExp” is dynamically changed for each sample.
resolution	uint32	500	Physical pitch (nm) between neighbor spots.
Dataset DataType: compound	DataType	Example (bin1)	Description
x	uint32	71032	x coordinate in bin N.
y	uint32	103180	y coordinate in bin N.
count	uint8/uint16/ uint32	1	MID count at (x, y) when bin size is N. Data type for “count” is consistent with “maxExp” in the “Attributes.”

geneExp/binN/gene:

Dataset “gene” is a 1D array which stores the gene names, the starting row indexes in dataset “expression”, and row counts.

Dataset	DataType: compound	Data Type	Example (bin1)	Description
gene		S32	b'Gm16045'	Gene name.
offset		uint32	21	The starting row index in dataset “expression” for the gene. In this example, the gene expression data for gene “Gm16045” starts from row 21 in the dataset “expression.”
count		uint32	2	Row count. In this example, expression data for gene “Gm16045” is recorded in row 21 and 22 (2 rows) in the dataset “expression.”

wholeExp/binN:

Dataset “binN” is a 2D array (matrix) which stores the MID count and gene type count at each spot.

Dataset Attributes	DataType	Example (bin1)	Description
number	uint64	22879557	Number of non-zero spots in the dense matrix.
minX	uint32	59820	Minimum x coordinate in bin N.
lenX	uint32	13221	Length of x.
minY	uint32	102086	Minimum Y coordinate in bin N.
lenY	uint32	18454	Length of y.
maxMID	uint8/uint16/ uint32	2155	Maximum MID count in a spot. Data type for “maxMID” is dynamically changed for each sample.
maxGene	uint16	846	Maximum gene type count in a spot.
resolution	uint32	500	Pitch (nm) between neighbor spots.

Dataset	DataType: 2D array (XxY), compound	Data Type	Example (bin1)	Description
MIDcount		uint8/uint16/ uint32	1	MID count in the spot. The spot coordinate can be identified from the row and column index of the 2D matrix plus the “minX” and “minY” specified in the attributes. Data type for “MIDcount” is consistent with “maxMID” in the “Attributes.”
genecount		uint16	1	Gene count in the spot. The spot coordinate can be identified from “Attributes” and the indexes of the 2D array.

stat/gene:

Dataset “gene” is a 1D array which stores the MID count and spatial pattern enrichment score (E10) of each gene. The array is order by the MID count in descending order.

Dataset Attributes	DataType	Example	Description
maxE10	float32	65.53	Maximum E10 score.
minE10	float32	0.	Minimum E10 score.
cutoff	float32	0.1	Threshold for filtering spots that will be used for computing E10. In this example, 0.1 means that the spots whose MID count is in the top 10% are used for calculating the spatial enrichment score.
Dataset Data Type: compound	DataType	Example	Description
gene	S32	b'Ptgds'	Gene name.
MIDcount	uint32	229502	MID counts for the gene.
E10	float32	65.53	The spatial pattern enrichment score (E10) for the gene.

The distinctions of each SAW output GEF files are explained in Table 2-3.

Table 2-3 SAW Output GEF Files Description

GEF Name	SAW Pipeline	Example	Description
SN.raw.gef	count	SS200000135TL_ D1.raw.gef	count output raw GEF, it only includes geneExp Group for the bin size of 1. The origin of expression matrix has been calibrated to (0,0).
SN.gef	tissueCut	SS200000135TL_ D1.gef	tissueCut output full GEF file. It contains geneExp Group and wholeExp Group for the bin size of 1, 10, 20, 50, 100, 200, and 500. SN.gef is also the only one that includes stat Group. The origin of expression matrix has been calibrated to (0,0), and its offsets are the same with SN.raw.gef. SN.gef is the input file for visualization.
SN.tissue.gef	tissueCut	SS200000135TL_ D1.tissue.gef	tissueCut output GEF file for the tissue-covered region. It only includes geneExp Group for the bin size of 1. The coordinates in the matrix and the offsets are all same with SN.raw.gef.

2.4 Gene Expression Matrix

Gene expression matrix stores genes spatial expression data. SAW generates multiple gene expression matrix files in the workflow, the basic format requires four columns with a header row that show the column names. The four columns are gene name, x coordinate, y coordinate, and MID count. The origin of **tissueCut** generated gene expression matrices have been calibrated to (0, 0). The header of expression matrix for maximum area enclosing rectangle region has six annotation rows start with “#” before the column rows. The header field names and field types are described in Table 2-4.

Table 2-4 Gene Expression Matrix Header Fields Description

Fields	Data Type	Example	Description
#FileFormat	string	GEMv0.1	Gene expression matrix file format version.
#SortedBy	string	None	Gene expression matrix sorting strategy. Valid values: “geneID”, “x”, “y”, “MIDCount”, “None”.
#BinSize	uint16	1	(Please check 1.3 Terminologies and Concepts Bin)
#StereoChip	string	SS200000135TL_D1	STOmics Gene Expression Chip SN.
#OffsetX	uint32	1	X coordinate of the origin before calibration.
#OffsetY	uint32	1	Y coordinate of the origin before calibration.
geneID	string	Cr2	Gene name.
x	uint32	16809	X coordinate of the spot.
y	uint32	8546	Y coordinate of the spot.
MIDCount	uint32	1	Number of MIDs at (x, y) for the gene in the corresponding row.

Example of GEM:


```

...
#FileFormat=GEMv0.1
#SortedBy=None
#BinSize=1
#StereoChip=SS200000135TL_D1.gem
#OffsetX=0
#OffsetY=0
geneID  x      y      MIDCount
Ccl27   4665   19736   1
CR974586.5  10207  15257   1
CR974586.5  4707   18336   1

```

2.5 Image Pyramid

The image pyramid model is a multi-resolution hierarchical model that is used to store and display images in different resolutions. For the same field of view, the layer of the image pyramid that is closest to the bottom includes the most detailed information and has the largest scale. **register** pipeline performs the down-sampling step on the registered image, and the resulted images are layered to construct a pyramid. For each resolution layer, the intact registered image is split into 256 pixels x 256 pixels tiles. If the size of a layer is smaller than 256 x 256, the image will then remain intact. The suffix of the file can be “.ssDNA.rpi” or “.conA.rpi” which indicates a cell nucleus or membrane stained image were used/registered respectively. Schematic diagram of image pyramid:

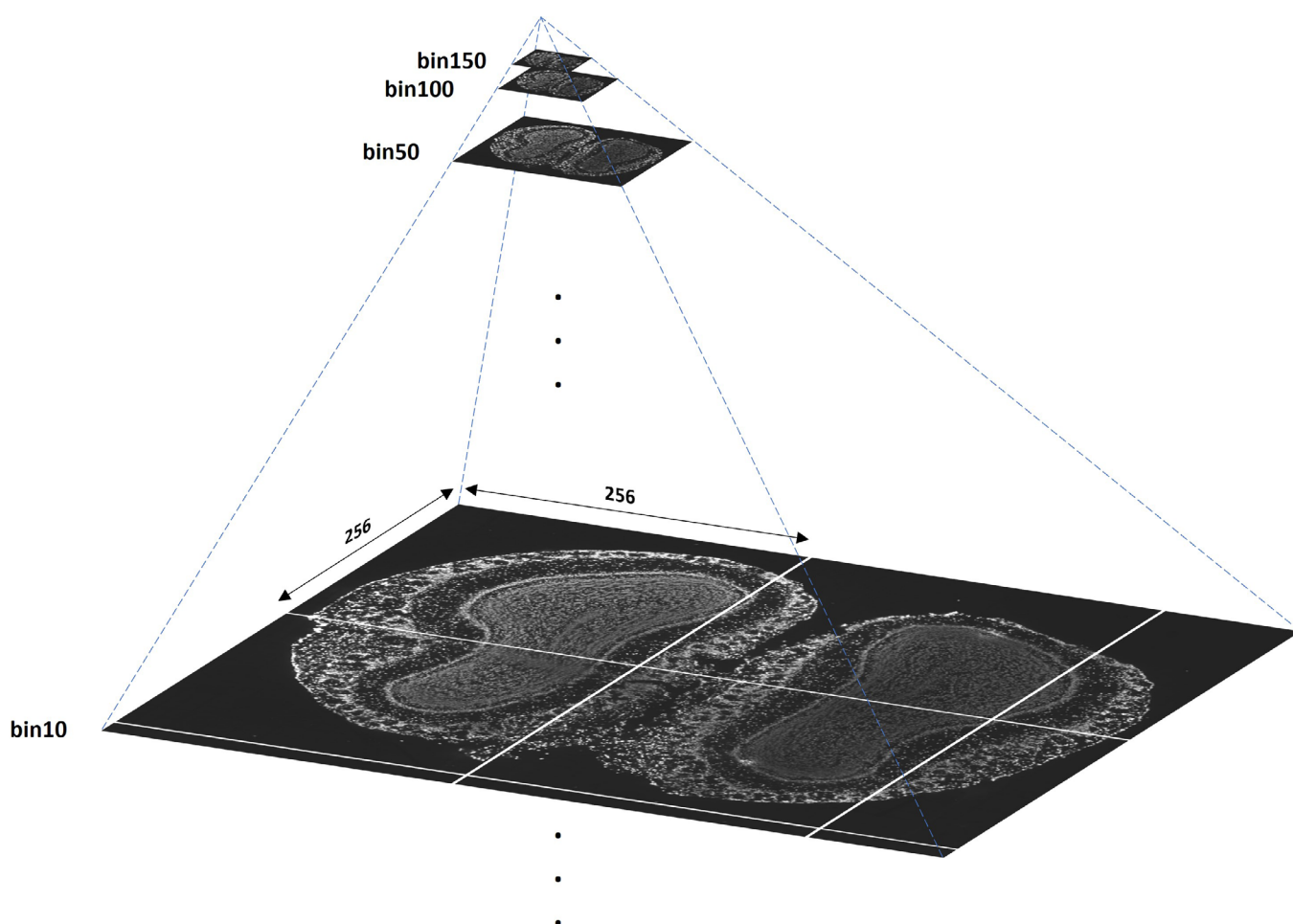

References

1. BGIResearch/SAW. Accessed October 13, 2021. <https://github.com/BGIResearch/SAW>
2. Sequence Alignment/Map Format Specification.; 2021. Accessed May 21, 2021. <https://github.com/samtools/hts-specs>.

Contact Us

BGI-Research, Shenzhen

stereomap.cngb.org/SAP

Email: support@stereomics.com